


Insufflator

Dyonics® Access 40 Insufflator

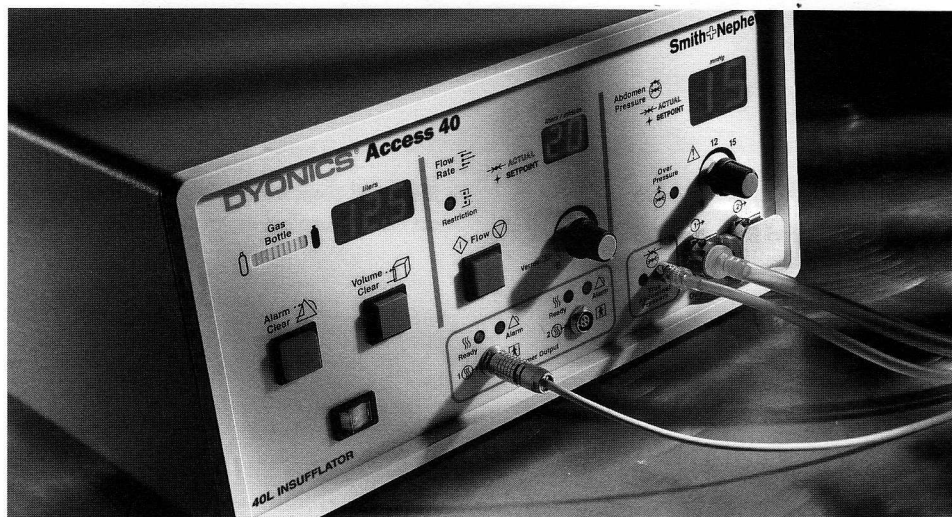

For fast, reliable insufflation, choose the Dyonics® Access 40 Insufflator which offers a unique method, using the TAP™ line, to monitor intra-abdominal pressure. With TAP™, Access 40 continuously calculates and measures abdominal pressure while the system is pumping, so you can rely on a faster response in the event of leaking during surgery.

Access 40's other high performance features include two CO₂ outputs that each pump 20 liters per minute, saving insufflation time. The system can be used with or without heat for flexibility in setting optimal conditions for the patient; heat reduces recovery time by lowering the risk of hypothermia during laparoscopy. An optional, reusable in-line gas warmer and disposable tubing make Access 40 more cost-effective than other systems.

- » Offers a unique method of monitoring intra-abdominal pressure using the TAP™ line while continuously pumping
- » Two CO₂ outputs for up to 40 liters per minute, providing optimum flow rates
- » Faster insufflation and reinsufflation in the event of loss of pneumoperitoneum or leaks
- » Reusable in-line warmer greatly reduces cost of consumables
- » Can be used with or without heat

REF	Name	Description	Price
Dyonics Access 40 Insufflator			
7209063	40L Insufflator	100 – 240 VAC, U.S. Kit	\$8,700
7206005	40L Insufflator	100 – 240 VAC, U.K. Kit	
7206006	40L Insufflator	100 – 240 VAC, French Kit	
7206007	40L Insufflator	100 – 240 VAC, German Kit	
7206008	40L Insufflator	100 – 240 VAC, Spanish Kit	
7206009	40L Insufflator	100 – 240 VAC, Italian Kit	
7205868	In-Line Gas Warmer	Reusable	895
7205861	Insufflation Tubing	Disposable (10 kits/box)	136
7205836	TAP™ Tubing Set	Disposable (10 kits/box)	110
4270	High Pressure Hose/Clamp	For connecting extra CO ₂ source	151
7208116	Disposable Pneumoperitoneum Insufflation Needle	14G x 12 cm length (10/box)	128
7205972	DISS Connector with 6 foot Hose		*

* Please contact your local authorized Smith & Nephew Representative for price information.

For Customer Service in the U.S. please call **1-800-343-5717**. Outside the U.S. please contact your local authorized Smith & Nephew Representative.